

CHARM SCIENCES, INC.

ROSA OCHRATOXIN QUANTITATIVE TEST

TABLE OF CONTENTS	PAGE
GENERAL INFORMATION	1
PREPARATION OF TESTING MATERIALS AND EQUIPMENT	2
EXTRACTION PROCEDURES	3
SAMPLE PREPARATION FOR QUANTIFICATION	3
TEST PROCEDURES	4
SUPPLEMENTAL ANALYSIS	6
REPORTING AND CERTIFYING TEST RESULTS	6
STORAGE CONDITIONS AND PRECAUTIONS	7
EQUIPMENT AND SUPPLIES	7
REVISION HISTORY	8

GENERAL INFORMATION

ROSA Ochratoxin Quantitative Test is an immunoreceptor assay utilizing ROSA (Rapid One Step Assay) lateral flow technology. Ochratoxin is extracted from the sample using 70% methanol in water. Ochratoxin A interacts with colored beads in the lateral flow test strip and the color intensity in the test and control zones is measured by the ROSA-M Reader and interpreted as parts per billion (ppb) ochratoxin A.

The instructions presented in this document cover only the procedure for performing the analytical test for official inspections. For questions regarding this procedure, contact Dr. Ajit Ghosh of the Technology and Science Division by phone at 816-891-0417 or email at Ajit.K.Ghosh@usda.gov.

Refer to the current policies and/or instructions issued by the Policies, Procedures, and Market Analysis Branch (PPMAB) of the Field Management Division for information on use of this test kit in official inspections including sampling, general sample preparation (e.g., grinding and dividing), reporting and certification of test results, laboratory safety, and hazardous waste management. For questions regarding these policies and/or instructions, contact Patrick McCluskey of PPMAB by phone at 816-659-8403 or email at Patrick.J.McCluskey@usda.gov.

Approved Test Kit Information

Test Kit Vendor:	<i>Charm Sciences, Inc.</i> 978-687-9200
Test Kit Name:	ROSA Ochratoxin Quantitative Test
Product Number:	LF-OCHRAQ-G
Effective Date of Instructions:	03/30/2015
Instructions Revision Number	0
Conformance Range:	5 – 100 ppb
Number of Analyses to Cover Conformance Range:	2
Type of Service:	Quantitative
Supplemental Analysis:	Yes
Approved Commodities:	Wheat, barley, corn, corn gluten meal, malted barley, oats, rye, sorghum, and soybean meal.
Extraction method:	Shake 50 gram sample with 100 milliliters (mL) 70% methanol/30% distilled or deionized water (v/v) for 1 minute.
Test Format:	Lateral flow strip
Detection Method:	ROSA-M Reader, Model LF-ROSAREADER-M-NB

PREPARATION OF TESTING MATERIALS AND EQUIPMENT

a. Test Strips:

Remove from the container only the number of test strips to be used in 1 day, document time of removal. Keep these test strips at room temperature during daily use for up to 12 hours and unused test strips should be discarded.

b. OCHRAQ Dilution Buffer:

- (1) Dispense buffer into a clean micro-centrifuge tube and label for each sample to be tested.
- (2) Use pre-dispensed buffer tubes and buffer solution at room temperature (18 to 30 °C).

c. Extraction Solvent [70% Methanol/30% Water (v/v)]:

The extraction solvent used in the method is a methanol/water mixture consisting of 70% methanol (reagent grade or better) and 30% distilled or deionized water (v/v).

- (1) Using a 1000 mL graduated cylinder, measure 700 mL methanol and place it into a clean carboy with spigot.
- (2) Using a 500 mL graduated cylinder, measure 300 mL distilled or deionized water and add to the methanol and shake until it is completely mixed.
- (3) Label the container stating the mixture 70% methanol/30% water (v/v), date of preparation, and initials of technician who prepared the solution.
- (4) Store this solution at room temperature in a tightly closed container until needed. Mix again before use.

NOTE: To prepare smaller or larger amounts of solution use the ratio of 7 parts methanol to 3 parts distilled or deionized water.

d. Negative Control:

Prepare negative control by adding 100 microliters (µL) extraction solvent to 1.0 mL OCHRAQ Dilution Buffer in a micro-centrifuge tube. Cap, mix and label.

e. Positive Control:

- (1) Reconstitute the dry positive control (provided with test kit) by adding 300 µL extraction solvent followed by 3.0 mL OCHRAQ Dilution Buffer. Mix well; allow to stand for 10 minutes at room temperature before use, and mix again just before use.
- (2) Store reconstituted positive control refrigerated (0 to 7 °C) for up to 1 week or aliquot (at least 1.5 mL) to clean micro-centrifuge tubes, label, and freeze within 6 hours of reconstitution at -15 °C or below for up to 2 months. Thaw slowly (overnight in refrigerator or with cool water) and shake well before use. Store thawed positive control refrigerated and use within 24 hours of thawing; DO NOT REFREEZE.

f. Reader and Test Strip Performance Testing:

- (1) Equipment Setup
ROSA-M Reader: Enter performance mode in ROSA-M Reader by selecting **OCHRA** channel in 3-line mode (**OCHRA** flashing) and sequentially pressing ESC, 5, ENTER. Follow ROSA-M Reader prompts to test calibration strips (LOWCAL and HIGHCAL) and controls (NEGCONTROL and POSCONTROL).
- (2) Test calibration strips daily to verify ROSA-M Reader performance. Calibration strips must test/perform in the specified ranges.
- (3) Test negative control and positive control weekly to verify test strip performance. Valid control ranges are:
 - (a) Negative Control: less than or equal to 1 ppb
 - (b) Positive Control: 4 to 10 ppb

If calibration strips or controls do not perform in specified ranges, discontinue use and contact Charm Sciences for assistance. Notify your monitoring field office or TSD with any documented information for quality control purposes.

g. ROSA Incubator:

ROSA Incubator must be clean and level. The ROSA Incubator temperature must be at 45 ± 1 °C (the temperature indicator should match the incubator temperature).

EXTRACTION PROCEDURES

- (1) Weigh 50 ± 0.2 grams ground samples into a clean extraction container.
- (2) Add 100 mL extraction solvent.
- (3) Shake vigorously for 1 minute (use within 30 minutes).
- (4) Allow sample to settle for 1 minute to obtain a clear sample extract.
- (5) Transfer 1 to 1.5 mL extract into a clean micro-centrifuge tube, label, and centrifuge for 10 seconds to obtain clarified extract (use within 2 hours).
- (6) Repeat for additional samples.

SAMPLE PREPARATION FOR QUANTIFICATION

This test kit uses different testing sensitivity ranges (Diluted Extract and Second Diluted Extract) for reporting ochratoxin measurements for grain and commodities.

a. Sample Preparation of filtered Diluted Extract for 5 to 30 ppb quantitation:

- (1) Pipet 1.0 mL OCHRAQ Dilution Buffer into a clean micro-centrifuge tube.

- (2) Pipet 100 μ L clarified extract to micro-centrifuge tube containing 1.0 mL OCHRAQ Dilution Buffer, cap, mix (5 times inverting up and down), and label. This tube contains the Diluted Extract.
- (3) Filter each Diluted Extract.
 - (a) Draw Diluted Extract into 1 mL syringe and pass through Minisart RC15 syringe filter.
 - (b) Collect the filtered Diluted Extract in a clean micro-centrifuge tube and label.
- (4) Repeat for additional samples.
- (5) Use filtered Diluted Extract (use within 6 hours after preparation) as your test sample in Sample Analysis found in Test Procedures section (page 4).

b. Sample Preparation of Second Diluted Extract for 20 to 100 ppb quantitation:


- (1) Pipet 1.0 mL OCHRAQ Dilution Buffer into a clean micro-centrifuge tube.
- (2) Pipet 300 μ L filtered Diluted Extract to micro-centrifuge tube containing 1.0 mL OCHRAQ Dilution Buffer, cap, mix (5 times inverting up and down), and label. This tube contains the Second Diluted Extract.
- (3) Use Second Diluted Extract (use within 6 hours after preparation) as your test sample in Sample Analysis found in Test Procedures section (below).

TEST PROCEDURES

a. Sample Analysis:

- (1) Check that the ROSA Incubator temperature is 45 ± 1 °C.
- (2) Label test strip(s) to identify sample.
- (3) Place test strip in the ROSA Incubator with the flat side facing upward.
- (4) Hold the test strip flat in the ROSA Incubator and use tab to expose sample compartment by peeling tape back to “Peel to Here” line.

Avoid lifting the test strip and sponge under tape and bending back the white wick and sponge under the tape.

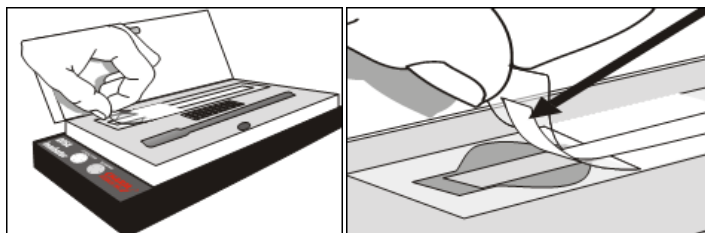

- (5) Hold the pipet vertically and slowly pipet 300 μ L test sample (filtered Diluted Extract, Second Diluted Extract, or control) into the sample compartment at the ROSA Incubator line.
- (6) Reseal the tape over the sample pad compartment.

NOTE: When performing multiple tests using a ROSA Incubator:

- (a) Peel, pipet, and reseal before starting next strip.
- (b) Complete all test strips within 1 minute.
- (7) Close lid on the ROSA Incubator.
- (8) Incubate test strip(s) for 10 minutes.
- (9) Remove strip from the ROSA Incubator.

Do not squeeze sample compartment. Hold test strip vertically with sample compartment in the down position until interpreted.


- (a) Wipe foreign matter (dust, etc.) from the test strip(s).
- (b) Inspect and read test strip within 2 minutes of incubation completion.

When running multiple test strips in the ROSA Incubator, remove one strip for visual inspection and interpretation at a time.

- (c) Lower ROSA Incubator lid; do not re-latch.

b. Visual Inspection:

- (1) The test strip is **INVALID** if any of the following are observed:
 - (a) C (Control) line is missing.
 - (b) T1, T2 (Test) or C line is smeared or uneven.
 - (c) T1, T2, or C line is obscured by diluted extract or control.
 - (d) Beads do not flow past T1, T2 or C lines.

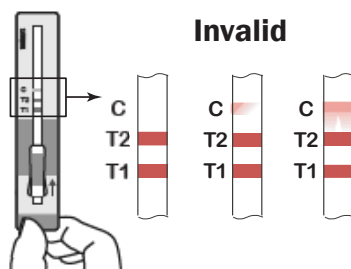

- (2) Do not put INVALID test strips in the ROSA-M Reader.
- (3) If test strip is INVALID, re-test the filtered Diluted Extract, Second Diluted Extract, or control.

c. Interpretation:

- (1) Insert a clean and valid test strip into the ROSA-M Reader. Slide the strip into the slot with the sample compartment in the up position until it stops.



- (2) Read results on **OCHRA** channel in 3-line mode (**OCHRA** flashing) using the appropriate **MATRIX**. If desired, enter Sample and/or Operator. Press ENTER to read.

- **MATRIX 00:** Assay of Diluted Extract for 5 to 30 ppb quantitation.
- **MATRIX 01:** Assay of Second Diluted Extract for 20 to 100 ppb quantitation.

For controls, see Reader and Test Strip Performance Testing in Preparation of Testing Materials and Equipment section (page 3).

- (3) **READING:** The number displayed is the concentration of ochratoxin (ppb) in the sample.

A “+” sign on a **READING** value indicates that the concentration of the sample is greater than the Sensitivity range. For example, a filtered Diluted Extract **READING** of “+30 ppb” indicates a value greater than 30 ppb. For quantitation of 20 to 100 ppb ochratoxin, prepare the Second Diluted Extract and use with another test strip.

A Second Diluted Extract **READING** less than 20 ppb indicates a value below the detection range. Re-test filtered Diluted Extract using another test strip for quantitation from 5 to 30 ppb ochratoxin.

A Second Diluted Extract **READING** greater than 100 ppb indicates that the concentration of the sample is greater than the sensitivity range of the sample dilution and is reported as “greater than 100 ppb”.

Note: Applicants may request qualitative certification in lieu of retesting of results outside of the filtered Diluted Extract or Second Diluted Extract test sample sensitivity ranges/concentrations.

SUPPLEMENTAL ANALYSIS

There are no supplemental analysis procedures approved for this test method. Sample results that report above the 100 ppb limit are certified as “greater than 100 ppb”.

REPORTING AND CERTIFYING TEST RESULTS

Refer to the current instructions issued by the Policies, Procedures, and Market Analysis Branch of the Field Management Division for reporting and certification of test results. For questions regarding these instructions, contact Patrick McCluskey (816-659-8403 or Patrick.J.McCluskey@udsa.gov).

STORAGE CONDITIONS AND PRECAUTIONS

a. Storage Conditions:

- (1) Store test strips refrigerated in tightly closed supplied container.
- (2) Store dilution buffer bottle and pre-dispensed micro-centrifuge tubes refrigerated.

Store reconstituted positive control refrigerated (0 to 7 °C) for up to 1 week or aliquot (at least 0.5 mL) to clean micro-centrifuge tubes, label, and freeze within 6 hours of reconstitution (-15 °C or below) for up to 2 months. Thaw slowly (overnight in refrigerator or with cool water) and shake well before use. Store thawed positive control refrigerated and use within 24 hours of thawing; DO NOT REFREEZE.

b. Precautions:

- (1) Test Strips
 - (a) To open test strip canister, remove and save plastic lid with foil lined foam insert to reseal container. Lift foil tab and peel foil seal off container. Discard foil seal.
 - (b) In high humidity, limit condensation by opening container after it has warmed to room temperature, estimated between 25 to 30 minutes from the time the container was removed from the refrigerator.
 - (c) Inspect/verify desiccant indicator. Beads inside desiccant packets should be blue. Do not use test strips if the blue beads have turned purple or pink
- (2) Use OCHRAQ Dilution Buffer supplied with each test kit only.
- (3) Do not use the test kits beyond the noted expiration date.
- (4) Debris on test strips may alter the reader optics. Keep equipment clean. Wipe dust and liquid off test strips before inserting into reader.
- (5) ROSA Incubator must be clean and level. ROSA Incubator temperature must be 45 ± 1 °C. The temperature indicator should match the ROSA Incubator temperature. A daily thermometer check is recommended. Keep ROSA Incubator lid lowered, but not latched unless performing test procedure. ROSA Incubator may take 10 minutes to reach proper temperature depending on ambient temperature.

EQUIPMENT AND SUPPLIES

a. Test Strips

- (1) LF-OCHRAQ-20K
 - (a) 1 container of 20 OCHRAQ test strips
 - (b) 1 OCHRAQ Grain Positive Control
 - (c) 1 OCHRAQ Dilution Buffer

- (2) LF-OCHRAQ-100K
 - (a) 1 container of 100 OCHRAQ test strips
 - (b) 1 OCHRAQ Grain Positive Control
 - (c) OCHRAQ Dilution Buffer
- (3) LF-OCHRAQ-500K
 - (a) 5 containers of 100 OCHRAQ test strips
 - (b) 5 OCHRAQ Grain Positive Controls
 - (c) 5 OCHRAQ Dilution Buffers

b. Materials required but not provided

- (1) 100 μ L pipet and pipet tips
- (2) 300 μ L pipet and pipet tips
- (3) 1000 μ L fixed volume pipet or 100 to 1000 μ L variable volume pipet and pipet tips
- (4) 100, 500, and 1000 mL graduated cylinders
- (5) Balance
- (6) Deionized or distilled water
- (7) Methanol (reagent grade or better)
- (8) Micro-centrifuge tubes
- (9) Mini-centrifuge
- (10) Minisart RC15 syringe filters (Sartorius Minisart RC 15, Part No. 17762)
- (11) ROSA-M Reader
- (12) Printer for ROSA-M Reader (optional)
- (13) ROSA Incubator
- (14) Sample extraction Whirl-pak bags or containers
- (15) Sample grinder
- (16) Storage bottle
- (17) Syringes
- (18) Transfer pipets (optional)

REVISION HISTORY

Revision 0 (03/30/2015)